WEHMER
ORTHODONTIC & DENTAL PRODUCTS

SÆBUR SUPPLY
a division of WEHMER
ORTHODONTIC & DENTAL LLC

Wehmer Dual Wheel Trimmer Maintenance Tips #8


Water Line Connection Kit # 3890

Thank -you for purchasing this convenience option for your Wehmer model trimmer. The #3890 Connection Kit is designed to facilitate hooking up your Wehmer model trimmer to the water supply . Each trimmer installation should include its own ***water supply shut-off valve*** *to which this connection kit attaches.*

Plumbing Installation Instructions

1. The water line connection kit includes two 1/8" NPT male threaded connectors. If the water supply valve to which you are connecting has different sized fittings, you will need to adapt the valve with the appropriate reducers and couplings.
2. With the water supply turned "OFF", connect one end of the waterline to the supply valve. *The brass fitting on either end of the plastic tubing is a two-piece fitting. Unscrew the part with threads and insert that into the water valve or adapter. Use pipe dope or Teflon thread sealing tape to assure a leak-free connection. Tighten the fitting securely with a 7/16" open end wrench. The brass piece with the black plastic insert should still be on the tubing. Slide it back far enough so you can slip the tubing into the other half of the fitting on the valve. Make sure it bottoms out inside the brass piece. While pushing the tubing into the fitting, screw the brass piece with the black plastic insert onto the mating half on the valve. Tighten with a 7/16" wrench.* DO NOT OVERTIGHTEN.
3. If your trimmer is fitted with an electric solenoid water valve, simply thread the remaining 1/8" NPT connector into the inlet side of the water valve. Again, use of pipe dope or Teflon thread sealing tape will assure a leak-free assembly. Follow instructions for water supply end.
4. If your trimmer is *not* fitted with an electric water valve, you must remove the threaded nipple (*part A in the illustration*) on the inlet side of your flow regulator. Then.... screw the connector from the water line into the flow regulator and secure per above instructions.

Remember, turn off your main water supply valve at the end of each work day!

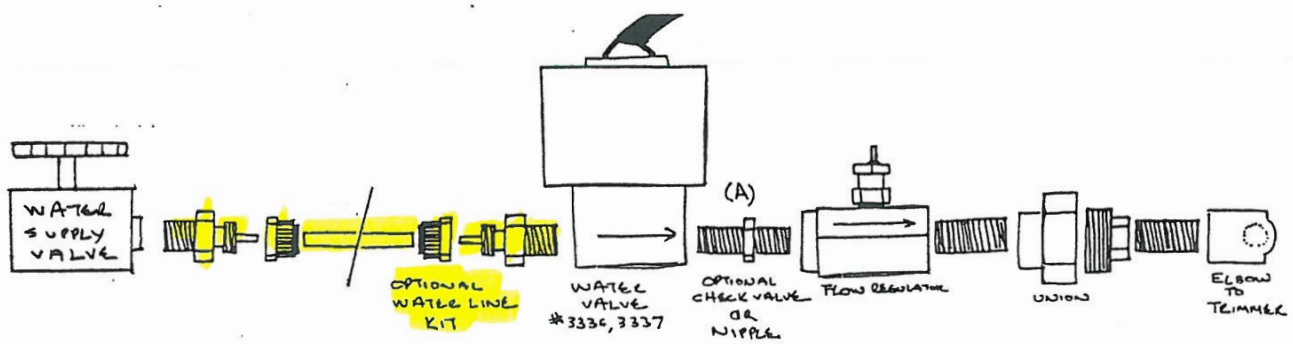

Questions? Wehmer Customer Service: 1-800-323-0229 \ 630-424-1877 \ customerservice@wehmer.com
FAX 630-424-1898 \ www.wehmer.com \ 483 W. Fullerton Avenue, Elmhurst, IL 60126